



**RMS**<sup>TM</sup>  
Home Infusion  
by RMS Medical Products

# 24 gauge needles...


The needle of choice to deliver viscous drugs up to 300ml/site/hr. Designed for sub-q with FREEDOM60<sup>®</sup> and Precision tubing, HlgH-Flo quality means less "ouch" and greater comfort throughout the infusion. When you require the maximum flow rate, select 24g HlgH-Flo. *Note: Standard 26 gauge HlgH-Flo is an excellent choice for viscous sub-q drugs up to 40ml/site/hr.*

needle tip has a backcut design which makes a tricuspid incision, resulting in less tearing of the skin

HlgH-Flo needle incision

other brand incision

Genuine 3M Tegaderm<sup>TM</sup> transparent dressing is included with all sets

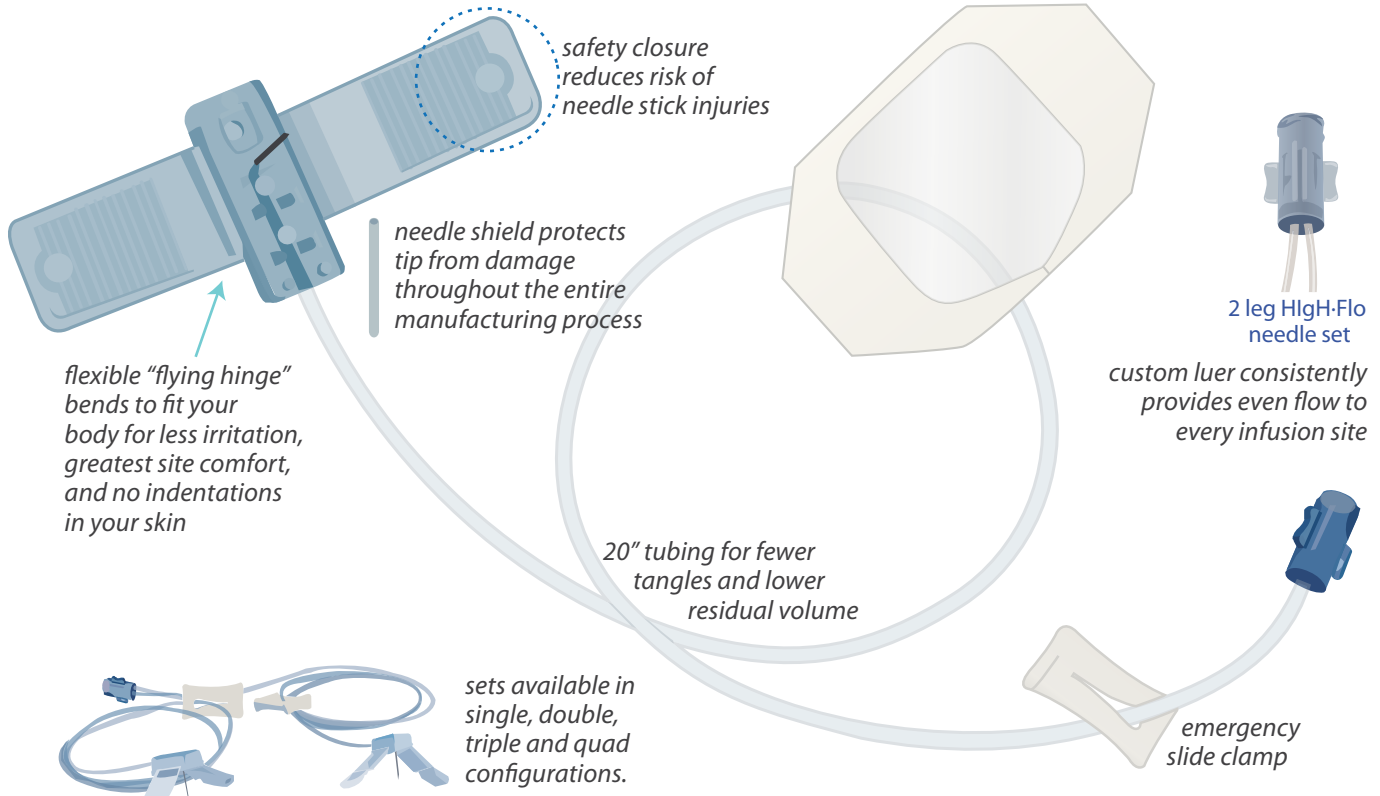

flexible "flying hinge" bends to fit your body for less irritation, greatest site comfort, and no indentations in your skin

needle shield protects tip from damage throughout the entire manufacturing process

safety closure reduces risk of needle stick injuries

20" tubing for fewer tangles and lower residual volume

sets available in single, double, triple and quad configurations.

2 leg HlgH-Flo needle set

custom luer consistently provides even flow to every infusion site

emergency slide clamp

## 24 Gauge Needle Set Ordering Information:

### 1-Needle Sets<sup>1</sup>

box of 20

Length	Item #
6mm	RMS12406
9mm	RMS12409
12mm	RMS12412
14mm	RMS12414

### 2-Needle Sets

box of 10

Length	Item #
6mm	RMS22406
9mm	RMS22409
12mm	RMS22412
14mm	RMS22414

### 3-Needle Sets

box of 10

Length	Item #
6mm	RMS32406
9mm	RMS32409
12mm	RMS32412
14mm	RMS32414

### 4-Needle Sets

box of 10

Length	Item #
6mm	RMS42406
9mm	RMS42409
12mm	RMS42412

### Related Options:



#### 24" Tubing Extender

Item #  
FEXT Box/20

#### Low Residual Volume Y-Connector

Item #  
LRVY Box/50



Y-Connector combines two needles, or allows infusing with two different needle lengths<sup>1</sup>.

<sup>1</sup> Note: To combine two RMS 1-Needle Sets, you must use the RMS LRVY. This is a part of the fluid path. Other y-connectors will affect medication delivery.



800 • 624 • 9600

rmsmedicalproducts.com

The FREEDOM60<sup>®</sup> Syringe Infusion System, FREEDOM60 Precision Flow Rate Tubing<sup>TM</sup> and HlgH-Flo Subcutaneous Safety Needle Sets<sup>TM</sup> are compliant with Medical Device Directive 93/42/EEC. RMS Medical Products is ISO 13485 certified.



RMS Medical Products  
Celebrating 35 Years  
Leadership & Service